
pure sense

CARE OF YOU

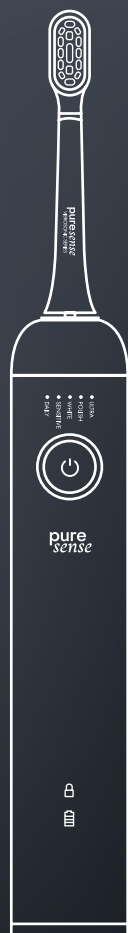Dental C

VIBROSONIC SERIES

PSVS300B • PSVS400B • PSVS500B • PSDV600

EN TABLE OF CONTENT

QUICK START 4-5

SAFETY PRECAUTIONS 6-7

USAGE..... 8

Parts..... 8

Before First Use 8

Charging 8

Brushing Instructions 8

Brushing Technique 8

Toothbrush Modes 9

Memory Function 9

Lock Function 9

MAINTENANCE10

Cleaning instructions10

UV sanitizer10

Storage10

SPECIFICATIONS11

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HURTIGSTART 4-5

SIKKERHET18-19

BRUK20

Deler 20

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Lading 20

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Børsteteknikk 20

Tannbørstемoduser 21

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Låsefunksjon..... 21

VEDLIKEHOLD22

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UV-desinfiseringsenhet22

Oppbevaring22

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SNABBSTART 4-5

SÄKERHET12-13

ANVÄNDNING14

Delar14

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Laddning14

Borstningsinstruktioner14

Borstningsteknik14

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Minnesfunktion15

Låsfunktion.....15



UNDERHÅLL.....16

Rengöringsinstruktioner.....16

UV-sterilisering16

Förvaring.....16

SPECIFIKATIONER 17

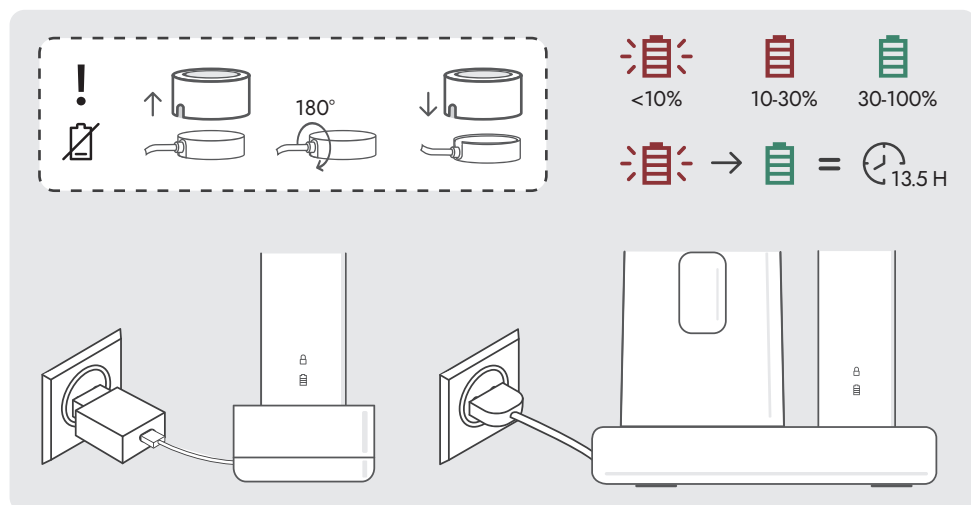
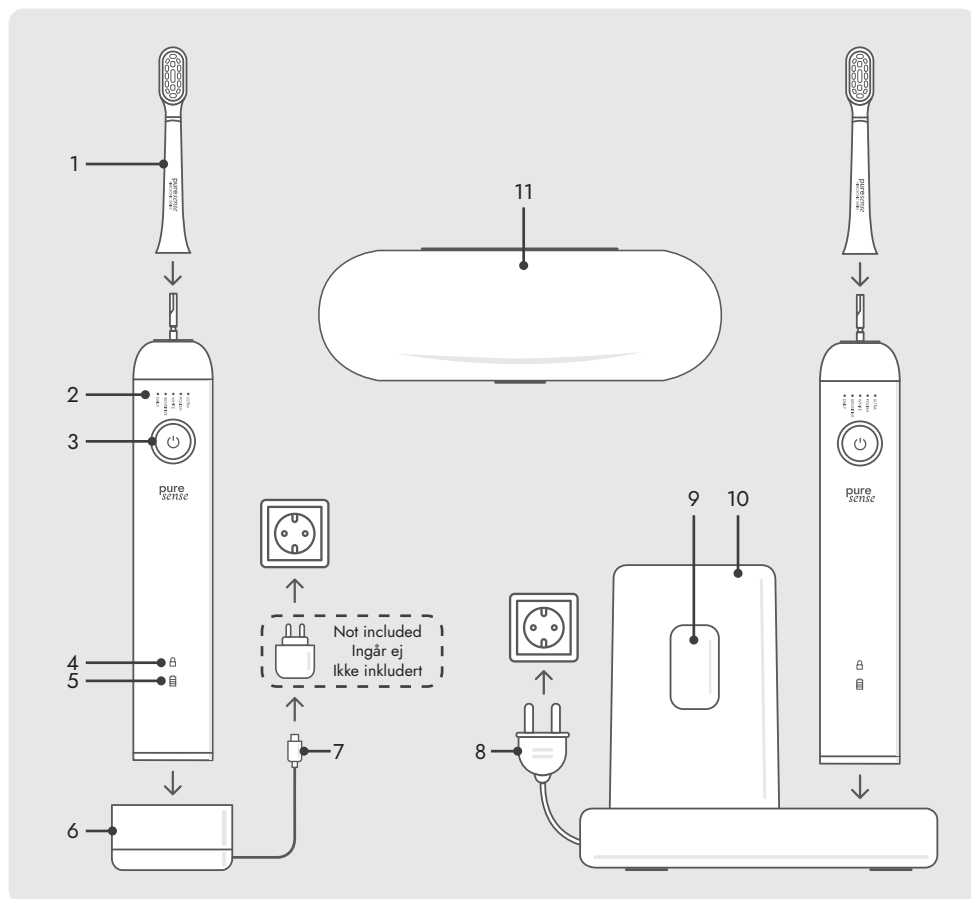
2 MIN

2 MIN

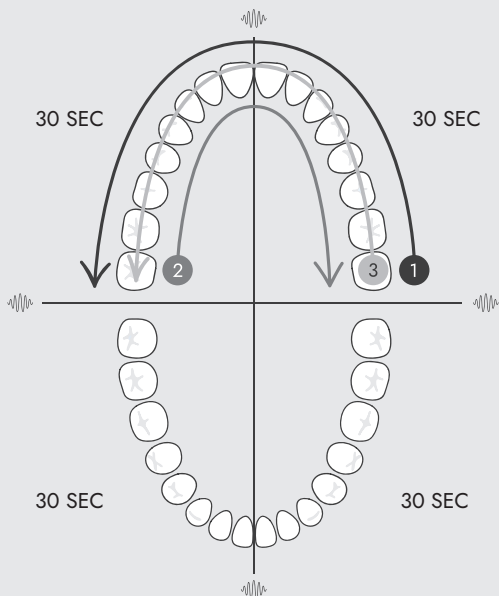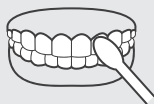DAILY


SENSITIVE

WHITE

POLISH

ULTRA

OFF

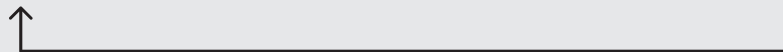Lock / Lås

1) Start the toothbrush

1) Starta tandborsten

1) Starte tannbørsten

2) Press again within 1 sec.
Hold for 5 sec.

2) Tryck igen inom 1 sek.
Håll in 5 sek.

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Hold inne i 5 sek.

SAFETY PRECAUTIONS

READ ALL OF THE INSTRUCTIONS BEFORE USE AND SAVE THEM FOR FUTURE REFERENCE

WARNING!

Risk of personal injury or product and/or property damage if you do not follow the precautions.

USAGE PRECAUTIONS

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Usage by children under age 3 is not recommended.
- Children shall not play with the appliance.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Any repairs to the product must be done by qualified personnel. Do not try to repair the product by yourself.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- This product is intended only for oral hygiene and for brushing the teeth and between the teeth. Any other use may result in damage to the device or in injury.
- For toothbrushing only. Do not chew or bite on brush.
- This product is intended only for domestic use and is not suitable for commercial use.
- Periodically check the entire product for damage. A damaged or non-functioning product should no longer be used.
- Prior to each use, ensure the brush head

fits well. Discontinue use of the toothbrush if brush head no longer fits properly. Never use without a brush head.

- Never insert any object into any opening of the appliance.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- Only use the charger provided with the appliance.
- Replace brush head every 3rd month or sooner if brush head becomes worn.
- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Do not use the product while bathing or showering.

BATTERY PRECAUTIONS

- Disposal of a battery into fire or a hot oven, or mechanically crushing or cutting of a battery can result in an explosion.
- Leaving a battery in an extremely high temperature can result in an explosion or leakage of flammable liquid or gas; and a battery subjected to extremely low air pressure can result in an explosion or leakage of flammable liquid or gas.
- The unit may get hot during use, this is normal.
- A new battery with 100% capacity reduces its capacity by 20% after approximately 300 charging cycles.
- Disconnect the appliance from charging as soon as it is fully charged.
- Charge the battery every 3rd month if you do not use the product for a longer period of time. This is to ensure the battery performance.
- This product contains a lithium battery. When a device with a lithium battery is damaged (either once or several times repeatedly), the lithium-ion cells can be affected by a so-called 'thermal rush'. Thermal rush means that energy is released in the form of heat inside the battery and fire occurs quickly. It becomes a self-generating process that takes

place explosively and is very difficult to extinguish.

- Do not apply excessive force on the product (do not hit or drop the product).
- Do not charge the product if it is very warm, to avoid charging a warm battery.
- Do not cover the product.
- Always inspect the product and stop using it if leakage, strange odors, color changes, deformation or any unusual occurrences are detected. If any leakage goes into your eyes, never rub them - clean your eyes with water and get medical help.
- Never disassembly or refit the unit. Repairs to the equipment shall only be performed by qualified personnel. Improper repairs may place the user at serious risk.

MAINTENANCE PRECAUTIONS

- Do not use solvents, scouring agents, hard brushes, metal tools, or sharp objects to clean the toothbrush.
- Solvents can harm your health and may corrode plastic parts.
- Abrasive cleaning agents and tools can scratch surfaces, reducing durability.
- Do not clean the charging unit under running water.
- Do not allow water to enter the interior of the housing.
- If water penetrates the device, have it inspected by a qualified electrician before using it again.
- This product shall not be thrown away or incinerated but disposed of in accordance with your local regulations concerning chemical wastes.

UV LIGHT PRECAUTIONS

UV light sanitizer is only available for some models. See packaging to confirm what your model includes.

- Never look directly at the UV light when activated.
- Keep fingers, skin, and eyes away from the UV source.
- Use only toothbrush heads designed for the sanitizing charger/case to ensure proper UV coverage and fit.

- Avoid placing non-compatible items (e.g., jewelry, non-brush items) inside.
- Ensure the UV light automatically shuts off when you open the cover/case. If not, unplug it and contact an authorized service center. Stop using the sanitizing charger/case immediately if you notice cracks, or malfunctions.

SYMBOL	EXPLANATION
	Read all of the instructions before use and save them for future reference.
	This manual does not provide any warranty or related declarations. All information in this manual including the product specification, pictures and literatures is just for reference. NetOnNet AB reserves the rights to update contents of this manual without further notice.
	This manual is the copyright of NetOnNet AB. No part of this manual may be revised, copied or transmitted in any way without written permission from NetOnNet AB.
	The CE mark indicates that this product has been manufactured according to current regulations for electromagnetic fields, which means that it is safe to use if used according to safety instructions. The product is manufactured according to current regulations within the EU legislation, limiting the use of hazardous substances in electrical and electronic equipment.
	This product and the accessories packed together are the applicable product to the WEEE directive except batteries. Please recycle any materials in accordance with your local recycling regulations. Batteries should never be thrown away or incinerated but disposed of in accordance with your local regulations concerning chemical wastes.
	This product's packaging materials are recyclable and can be reused. Please recycle any materials in accordance with your local recycling regulations.
IPX7	Waterproof
	This symbol identify equipment meeting the safety requirements specified for Class II equipment according to IEC 61140.
	The exclamation point within an equilateral triangle is intended to alert the user to the presence of important instructions. If these instructions are not followed, there is risk of personal injury or product and/or property damage.

USAGE

PARTS

- 1. Brush head
- 2. Mode indicators
- 3. Power button
- 4. Lock indicator
- 5. Battery indicator
- 6. Charger*
- 7. USB-C connector
- 8. Power plug
- 9. UV lamp (on the inside)
- 10. Charger with UV sanitizer*
- 11. Travel case

* Some models include a UV sanitizer integrated with the charging station (10). Not all toothbrushes come with this feature, some has a charging unit without sanitizer (6). Please check your product packaging to confirm what your model includes.

BEFORE FIRST USE

- Carefully remove the toothbrush and its accessories from the packaging. Check that they are complete and undamaged.
- If you have any existing oral health conditions or are undergoing treatment, consult your dentist before use.
- Charge the toothbrush for 24 hours before first use.
- The vibrosonic technology is different from traditional electrical toothbrushes. It may take some time to get used to the technique. It is therefore recommended to start on SENSITIVE mode as a beginner.

CHARGING

- When the battery indicator (🔋) flashes in red, the battery level is less than 10% and the toothbrush must be charged.
- 1. Plug the charger into an electrical outlet.
- 2. Place the toothbrush on the charger.
- 3. While charging, the indicator light flashes green. Once fully charged, the light will be solid green (🟢).
- 4. Unplug the charger after charging.

LIGHT COLOR	BATTERY LEVEL
Green light	30-100%
Red light	10-30%
Flashing red light	<10%

Note: If the toothbrush does not charge when placed on the charger, the charging plate may be installed upside down. Flip the charging plate and try again.

BRUSHING INSTRUCTIONS

- 1. Attach the brush head to the handle (there will be a short 0.5 mm gap).
- 2. When picking up the handle, it will automatically center the brush head.
- 3. Wet the bristles slightly and apply toothpaste.
- 4. Place the toothbrush in your mouth before turning it on to avoid splattering.
- 5. Turn it on and select a mode by pressing the power button (🔘) multiple times.
Note: To switch mode, you need to press the button **again within 5 seconds**. After 5 seconds, the toothbrush will turn off when pressing the power button.
- 6. Brush for a total of two minutes, spending 30 seconds per quadrant (the toothbrush will alert you when to switch).
- 7. The toothbrush will automatically turn off after 2 minutes, or you can press the power button to turn it off manually.
- 8. Spit out the toothpaste but avoid rinsing the mouth immediately (this helps retain fluoride for added protection).
- 9. Clean your toothbrush after each use by following the cleaning instructions.

BRUSHING TECHNIQUE

- Hold the brush at a 45-degree angle to your gums.
- Gently move the brush head tooth by tooth, letting it rest on each for a few seconds.
- Let the toothbrush do the work - avoid scrubbing back and forth like a manual brush.
- Use light pressure - pressing too hard can wear down enamel and irritate gums.
- Focus on all surfaces: Outer, inner, and chewing surfaces.

TOOTHBRUSH MODES

Press the power button (⏻) to cycle through modes. The active mode indicator will light up.

Note: To switch mode, you need to press the button again **within 5 seconds**. After 5 seconds, the toothbrush will turn off when pressing the power button.

MODE	FUNCTION
DAILY	For regular comprehensive cleaning. Gently removes plaque and food debris while maintaining enamel health. Ideal for routine morning/night brushing with balanced vibration intensity that cleans effectively without excessive abrasion.
SENSITIVE	Sensitive teeth protection. Reduced vibration intensity minimizes gum irritation, perfect for users with gum recession, tooth sensitivity, or post-dental procedures. Ensures thorough cleaning with extra comfort.
WHITE	Surface stain removal. Combines slightly stronger oscillations to eliminate coffee/wine/tobacco stains. Enhances natural whiteness through micro-polishing action while protecting enamel.
POLISH	Professional-grade polishing. Uses high-frequency pulsations to smooth enamel surfaces, reducing plaque adhesion.
ULTRA	Intensive deep cleaning. Maximum power setting with multi-directional motions to disrupt stubborn plaque. Recommended for occasional use (1-2 times per week), before special events or after consuming staining foods.

MEMORY FUNCTION

After a mode is used for more than 5 seconds, the toothbrush will memorize the selected mode. The next time it is turned on, it will automatically select the memorized mode to start.

LOCK FUNCTION

The toothbrush has a lock function to prevent accidental activation.

Activating the Lock:

1. Press the power button to start the toothbrush.
2. Release the button, then press and hold it again within 1 second for approximately 5 seconds.
3. The lock indicator will light up, confirming activation.
4. The indicator will turn off after a few seconds but will light up briefly if the power button is pressed while locked.

Deactivating the Lock:

1. Press and hold the power button for approximately 5 seconds until the indicator stops flashing.
2. When the indicator turns off, the lock is deactivated, and the toothbrush can be used normally.

MAINTENANCE

CLEANING INSTRUCTIONS

To maintain hygiene and ensure proper functioning, clean the toothbrush and brush head after each use.

Cleaning the Brush Head & Handle:

1. Remove the brush head from the handle.
2. Rinse both parts separately under running water.
3. Wipe them dry.

Cleaning the Charging Unit:

1. Unplug the charger from the electrical outlet before cleaning.
2. Wipe the exterior with a dry, soft cloth.
3. If necessary, use a slightly damp cloth to clean the housing, then dry it thoroughly before plugging it back in.

Cleaning the Travel Case:

1. Wipe the inside of the case with a dry or slightly damp cloth to remove moisture and toothpaste residue.
2. For deeper cleaning, use warm water, mild detergent, and a soft brush or cloth to gently scrub away buildup.
3. Leave the case open to air dry completely before storing your toothbrush again.

UV SANITIZER

Some models include a UV sanitizer integrated with the charging station. Not all toothbrushes come with this feature. Check your product packaging to confirm if your model includes a UV sanitizer.

How the UV Sanitizer Works:

The sanitizer effectively eliminates 99.9% of bacteria on the brush head(s) between uses. To maintain hygiene, it is recommended to sanitize brush heads before and after each use.

Using the Sanitizing Function:

1. Connect the charging/sanitizing unit into an electrical outlet.
2. Insert the brush head(s) to the small shafts on the charging/sanitizing unit. Ensure the bristles are facing the UV lamp.
3. Close the cover to start the sanitization process. The UV lamp will turn on

automatically and run for approximately 7 minutes.

4. The UV lamp will turn off automatically when the cycle is complete.
5. If the cover is opened during operation, the UV lamp will shut off immediately for safety.

Note: The sanitizing function activates automatically each time the cover is closed, when the appliance is plugged in.

Cleaning the Charging/Sanitizing Unit:

1. Unplug the charger from the electrical outlet before cleaning.
2. Wipe the exterior and interior with a dry, soft cloth.
3. If necessary, use a slightly damp cloth to clean the housing, then dry it thoroughly before plugging it back in.

STORAGE


- The toothbrush should be stored away from dust, dirt and water sources such as tubs and sinks, or in a travel case.
- The charger should be stored away from dust, dirt and water sources such as tubs and sinks.
- Brush head can be stored on toothbrush or separately put inside a travel case or other suitable location. Keep it away from dust and dirt.

SPECIFICATIONS


TOOTHBRUSH SPECIFICATIONS	
Model name	PSVS260B / PSVS260W
Battery type	Li-ion, 2000mAh
Battery life	180 minutes
Charging time	11.5 hours
Oscillation speed	180 times/minute
Vibration speed	Up to 31 000 strokes/min
IP class	IPX7
Handle size	31 x 29 x 183 mm
Handle weight	106 g

TRAVEL CASE SPECIFICATIONS	
Model name	PSTC50CB
Size	214 x 74 x 41 mm
Weight	106 g

CHARGER* SPECIFICATIONS	
Model name	PSVS260B
Input	5V  1A
Power	3W
IP class	IPX7
Cable length	100 cm
Size	48 x 48 x 26 mm
Weight	50 g

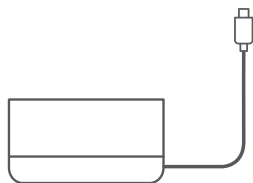CHARGER WITH UV SANITIZER* SPECIFICATIONS	
Model name	PSVS260BUV
Input	100-240V AC, 50-60Hz
Power	5.0W
Standby mode power consumption	0.3W
Auto enter standby mode period	5s
Off mode power consumption	No
Auto enter off mode period	No
IP class	IPX7
Sanitizing time	7 minutes
UV lamp wavelength	UVC: 270-280nm UVA: 395-405nm
UV lamp intensity	800μW/cm ²
UV lamp lifespan	> 1 300 hours
Cable length	135 cm
Size	67 x 127 x 114 mm
Weight	270 g

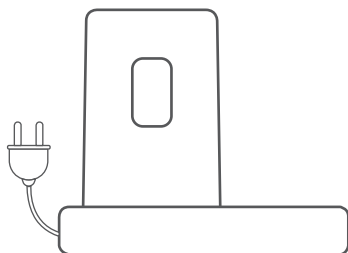* Some bundles includes a regular charger and some includes a charger with UV sanitizer. Check your product packaging to confirm which charger you have.

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Manufactured by NetOnNet AB