
*pure*sense

CARE OF YOU

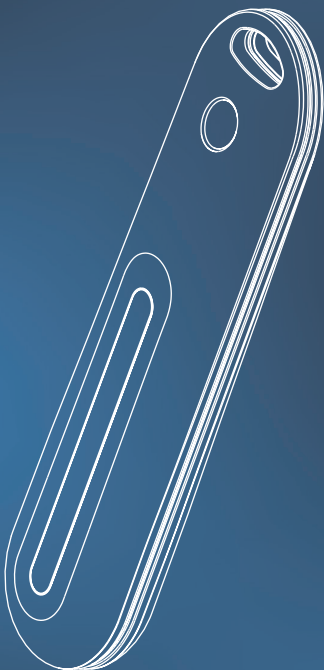Dental C

SONIC SERIES

PSTC500B

TABLE OF CONTENT

Safety precautions	4
Overview	6
Operation.....	7
Specifications.....	8

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SAFETY PRECAUTIONS

READ ALL OF THE INSTRUCTIONS BEFORE USE AND SAVE THEM FOR FUTURE REFERENCE.

WARNING! Risk of personal injury or product and/or property damage if you do not follow the precautions.

USAGE PRECAUTIONS


- Do not place or store the charger case where it can fall or be pulled into a bathtub, sink or toilet.
- Do not place or drop in water or other liquid. Do not reach for the charger case that has fallen into water or other liquid. Unplug immediately.
- Do not use while bathing.
- Unplug charger case before cleaning and assure that it is dry before plugging it into an electrical outlet.
- Never look directly at the UV light when activated.
- Keep fingers, skin, and eyes away from the UV source.
- Use only toothbrush heads designed for the sanitizing case to ensure proper UV coverage and fit.
- Avoid placing non-compatible items (e.g. jewelry, non-brush items) inside.
- Ensure the UV light automatically shuts off when you open the lid. If not, unplug it and contact an authorized service center. Stop using the sanitizing case immediately if you notice cracks, or malfunctions.

SYMBOL	EXPLANATION
	Read all of the instructions before use and save them for future reference.
	This manual does not provide any warranty or related declarations. All information in this manual including the product specification, pictures and literatures is just for reference. NetOnNet AB reserves the rights to update contents of this manual without further notice.
	This manual is the copyright of NetOnNet AB. No part of this manual may be revised, copied or transmitted in any way without written permission from NetOnNet AB.
	The CE mark indicates that this product has been manufactured according to current regulations for electromagnetic fields, which means that it is safe to use if used according to safety instructions. The product is manufactured according to current regulations within the EU legislation, limiting the use of hazardous substances in electrical and electronic equipment.
	This product and the accessories packed together are the applicable product to the WEEE directive except batteries. Please recycle any materials in accordance with your local recycling regulations. Batteries should never be thrown away or incinerated but disposed of in accordance with your local regulations concerning chemical wastes.
	This product's packaging materials are recyclable and can be reused. Please recycle any materials in accordance with your local recycling regulations.
	The exclamation point within an equilateral triangle is intended to alert the user to the presence of important instructions. If these instructions are not followed, there is risk of personal injury or product and/or property damage.

OVERVIEW

PARTS




1. Toothbrush charging
2. UV sanitizer for brush heads
3. Storage compartment
4. USB-C charging port
5. Power indicator

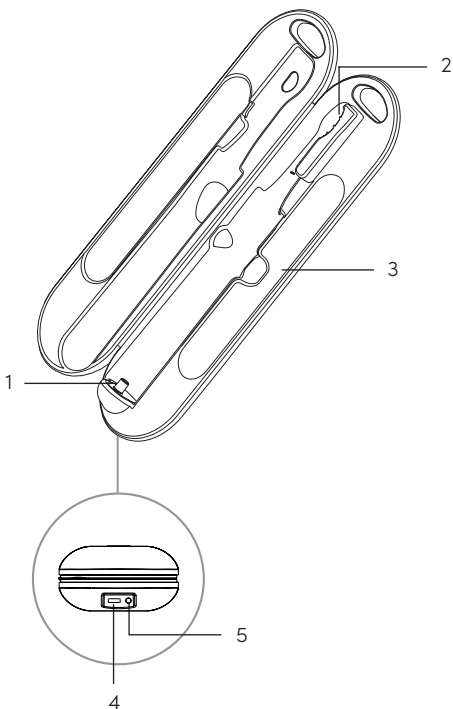OPERATION

USE THE CASE

1. Open the lid, place the suitable toothbrush (Sonic series toothbrush) into the charging travel case.
2. Plug a USB-C cable into the charging port of travel case, then plug another head into the adaptor or other USB power sources.
3. The power indicator will light in red. The toothbrush charge indicator on the handle will light up.

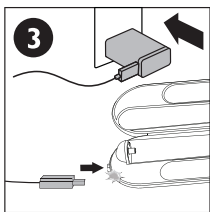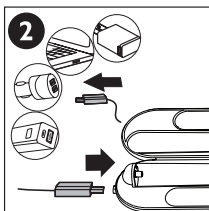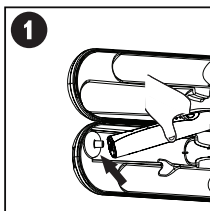4. Close the lid, UV sanitizer will start working automatically for around 7 minutes.
5. The power indicator lights up red, even if the toothbrush is fully charged.

MAINTENANCE

- Wipe the device with a dry cloth.
- The product should be stored away from dust, dirt and water sources such as tubs and sinks.
- Store the product out of reach of children.

SPECIFICATIONS

Model name	PSTC500B
Power supply	5V DC, 600mA
Sanitizing time	7 minutes
IP class	IPX4
Storage compartment	143 x 27 mm
Size	299 x 74 x 41 mm
Weight	211 g
Included in the package	- Travel case - USB-A to USB-C charging cable - User manual

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